


# FLEBOGRIF<sup>®</sup>

SET FOR VENOUS INSUFFICIENCY TREATMENT

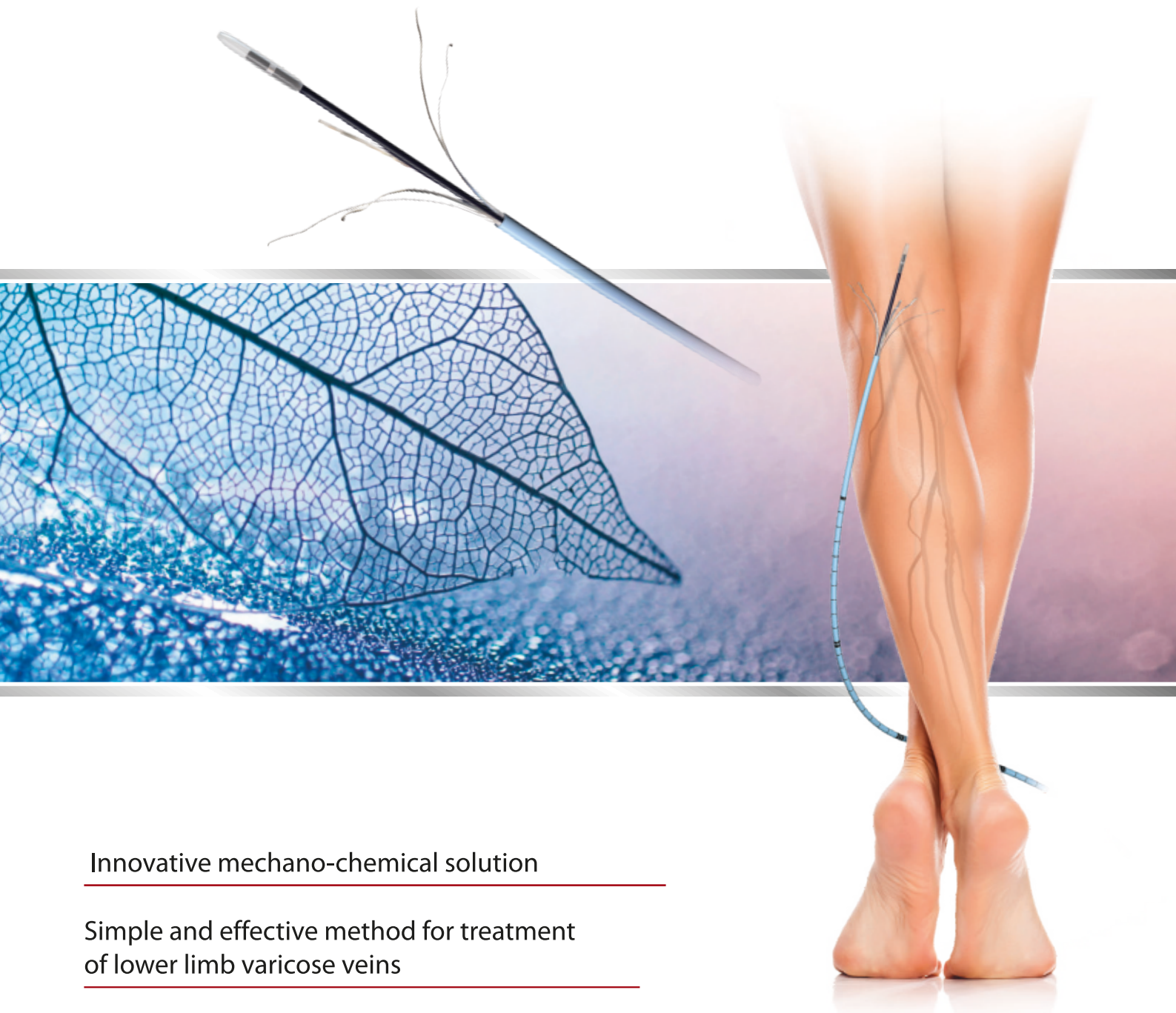

Innovative mechano-chemical solution

Simple and effective method for treatment  
of lower limb varicose veins

Clinically proven effectiveness

**BALTON**<sup>®</sup>



CE<sub>0197</sub>

# FLEBOGRIF®

SET FOR VENOUS INSUFFICIENCY TREATMENT

The set for varicose veins treatment Flebogrif is equipped with retractable scratching elements that combine the mechanical and chemical obliteration of the superficial venous trunks. The method consists in a simultaneous mechanical damage of the inner layer of the vein wall and chemical controlled administration of obliterating agent.


## Advantages:

- painless outpatient procedure
- no thermal action
- easy to use
- no other medical equipments (except USG) are needed
- quick procedure
- self-expandable scratching elements that adapt to the caliber of the vein
- excellent and long-lasting cosmetic effect
- quick recovery
- no tumescent anaesthesia

## Set elements:

- straight needle 18G
- guide wire J .035"
- vascular sheath with dilatator 6F
- catheter "Flebogrif" 6F
- three-way stopcock
- syringe 5 ml Luer-Lock (2 pcs.)

The effectiveness confirmed  
is 92% at 24-months  
observation\*

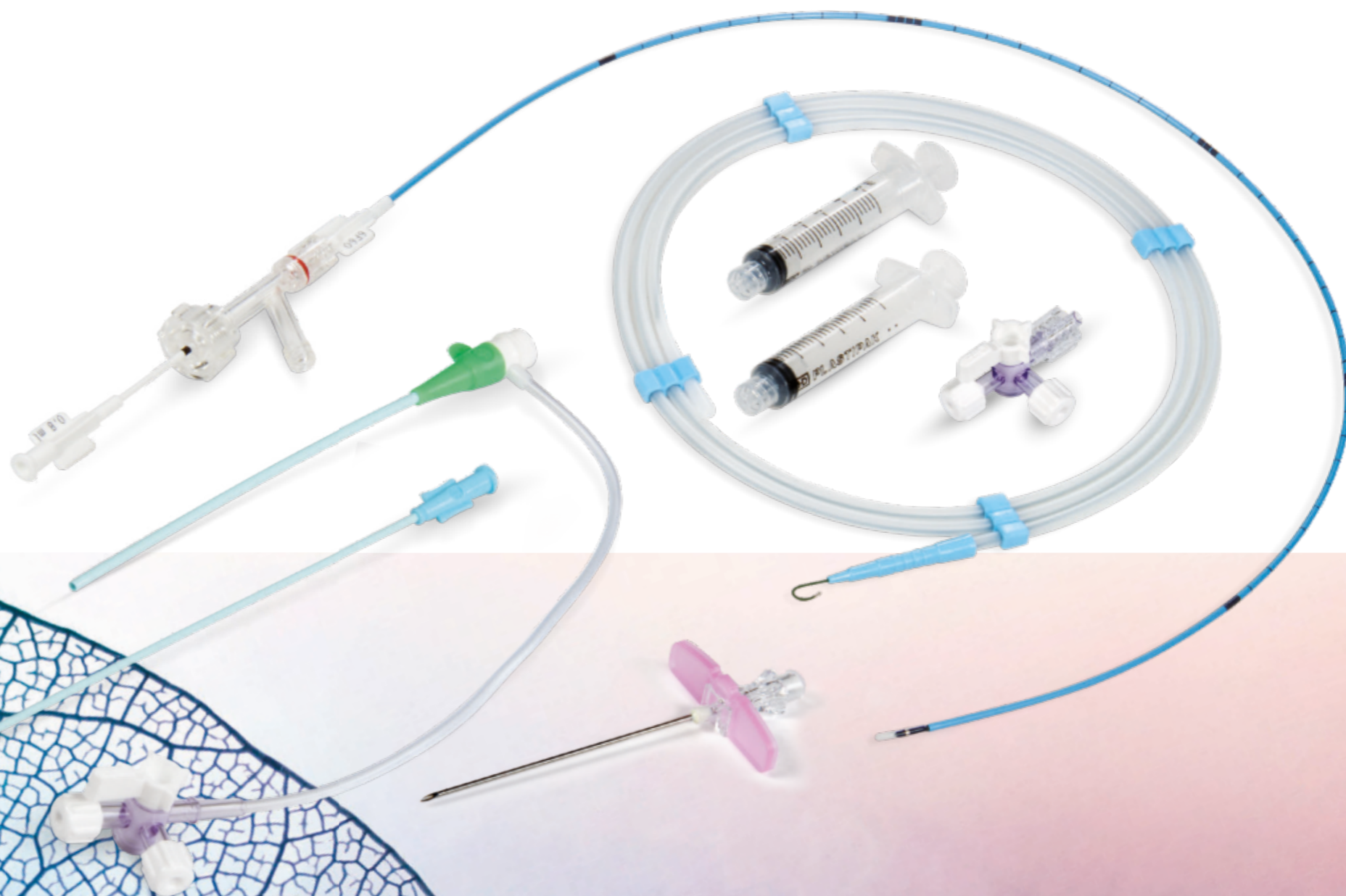

Catalogue No.	Catheter size	Catheter length
FLE6F90	6F	90 cm
FLE6F60	6F	60 cm

\*The novel minimally invasive mechano-chemical technique of the saphenous vein ablation. Our center experience: results of 24 months follow-up  
Marek Ilzecki, Piotr Terlecki, Stanislaw Przywara, Joanna Ilzecka, Shawn Dave, Tomasz Zubilewicz  
Acta Angiologica 2019; Vol. 25, No. 3, pp. 127-132 Doi: 10.5603/AA.2019.0008

# FLEBOGRIF®

SET FOR VENOUS INSUFFICIENCY TREATMENT



**Caution:**

The content of this document, including provided data, information, trademarks and logos are the sole property of BALTON. Any reproduction, whether in whole or in part, is strictly prohibited and will be considered a violation of BALTON's copyrights and other intellectual proprietary rights. ©2021 Balton Sp. z o.o. All rights reserved. This document and any images contained therein is intended solely for the use by healthcare professionals and BALTON's distributors.

This publication is for information purposes only and should not be construed as a commercial offer in the meaning of art. 66(1) of the Polish Civil Code.

All medical devices made commercially available by BALTON must be used in accordance with the enclosed instructions for use.

This publication presents only a small sample of our complete offering of medical devices. For further information please contact your Balton representative or email us at [balton@balton.pl](mailto:balton@balton.pl).



**Bioline Comercial Ltda.**

R. José Getúlio, 579 14º andar cj 141/146

Liberdade | SP | 01509 001

+55 11 3274 7600

[www.bioline.com.br](http://www.bioline.com.br)