


INTRODUCERS

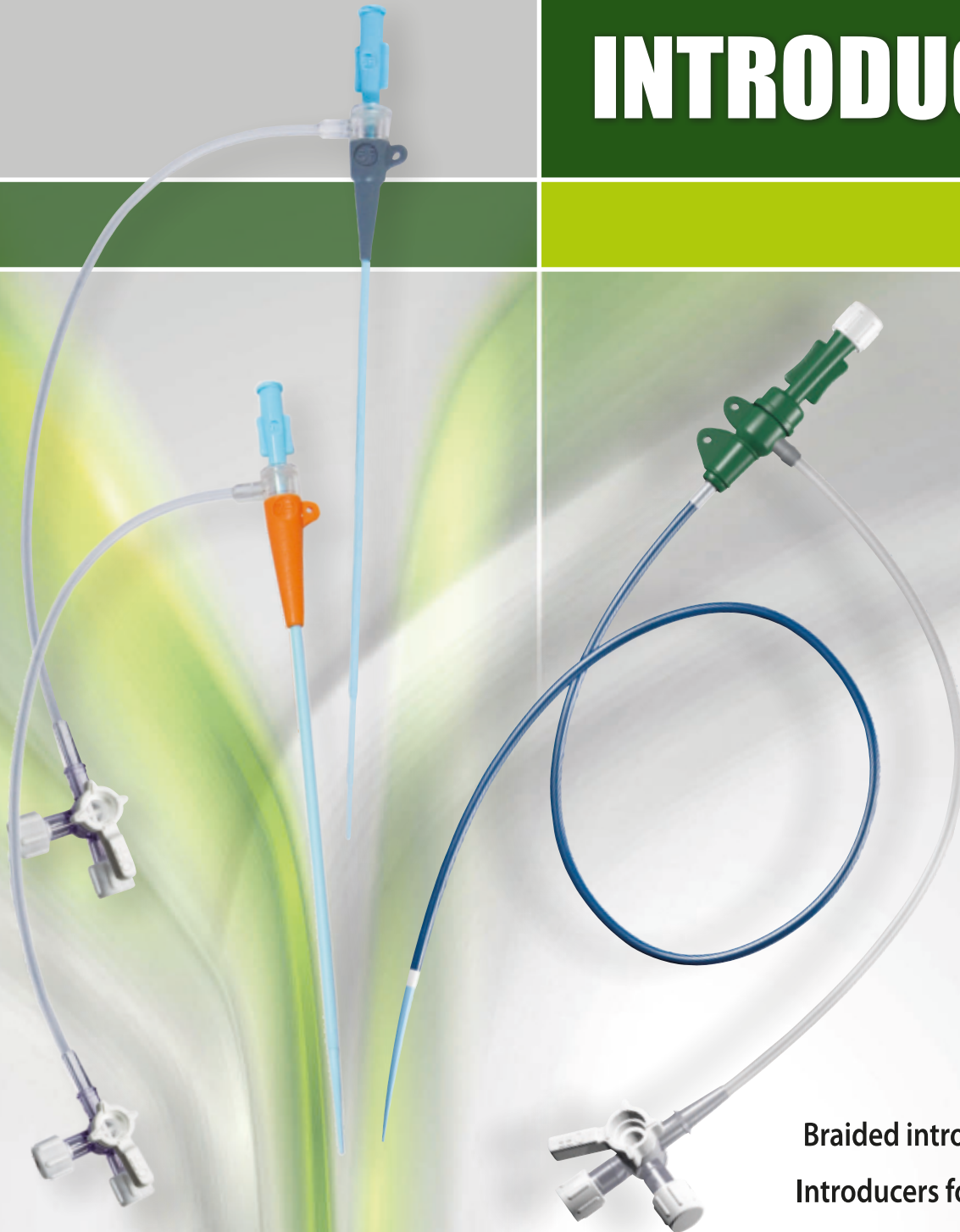

Braided introducers (Coil) ●

Introducers for cardiology ●

Introducers for radiology ●

Radial artery introducers for cardiology ●

BALTON
since
1980

BALTON®

CE 0197

INTRODUCERS

INTRODUCER

Set with hemostasis valve to introduce and exchange catheters and endocavitary electrodes

Set contains:

- sheath with a valve
- angiographic needle
- guide wire
- dilatator

Catalog Number	Sheath size / length	Cap colour	Dilatator size	Angiographic needle size	Guide wire size / length / type
INT3F	3F / 60 mm	white	3F	21G	.018" / 35 cm / STRAIGHT
INT4F	4F / 70 mm	red	4F	19G	.022" / 40 cm / "J"
INT5F	5F / 110 mm	grey	5F	18G	.035" / 40 cm / "J"
INT6F	6F / 110 mm	green	6F	18G	.035" / 40 cm / "J"
INT7F	7F / 110 mm	orange	7F	18G	.035" / 40 cm / "J"
INT8F	8F / 110 mm	blue	8F	18G	.035" / 40 cm / "J"
INT9F	9F / 110 mm	black	9F	18G	.035" / 40 cm / "J"
INT10F	10F / 110 mm	white	10F	18G	.035" / 40 cm / "J"
INT12F	12F / 120 mm	blue	12F	18G	.035" / 40 cm / "J"
INT14F	14F / 120 mm	blue	14F	18G	.035" / 40 cm / "J"
INT16F	16F / 120 mm	blue	16F	18G	.035" / 40 cm / "J"

Special order can contain introducer with radiopaque marker.
Model reference number: **INT8FM** – introducer with a marker.

INTRODUCER FOR CARDIOLOGY

Set with hemostasis valve to introduce and exchange catheters and endocavitary electrodes

Set contains:

- sheath with a valve
- angiographic needle
- guide wire
- dilatator

Catalog Number	Sheath size / length	Cap colour	Dilatator size	Angiographic needle size	Guide wire size / length / type
INT3FK	3F / 60 mm	white	3F	21G	.018" / 35 cm / STRAIGHT
INT4FK	4F / 70 mm	red	4F	19G	.022" / 40 cm / "J"
INT5FK	5F / 110 mm	grey	5F	18G	.035" / 40 cm / "J"
INT5FK23	5F / 230 mm	grey	5F	18G	.035" / 60 cm / "J"
INT6FK	6F / 110 mm	green	6F	18G	.035" / 40 cm / "J"
INT6FK23	6F / 230 mm	green	6F	18G	.035" / 60 cm / "J"
INT7FK	7F / 110 mm	orange	7F	18G	.035" / 40 cm / "J"
INT7FK23	7F / 230 mm	orange	7F	18G	.035" / 60 cm / "J"
INT8FK	8F / 110 mm	blue	8F	18G	.035" / 40 cm / "J"
INT8FK23	8F / 230 mm	blue	8F	18G	.035" / 60 cm / "J"
INT9FK	9F / 110 mm	black	9F	18G	.035" / 40 cm / "J"
INT10FK	10F / 110 mm	white	10F	18G	.035" / 40 cm / "J"

Special order can contain introducer with radiopaque marker.
Model reference number: **INT8FKM** – introducer with a marker.

LONG INTRODUCER

Set with hemostasis valve to introduce and exchange catheters and endocavitary electrodes

Set contains:

- sheath with a valve
- angiographic needle
- dilatator

Catalog Number	Sheath size / length	Cap colour	Dilatator size	Angiographic needle size
INT5FK50	5F / 500 mm	grey	5F	18G
INT5FK55	5F / 550 mm	grey	5F	18G
INT5FK60	5F / 600 mm	grey	5F	18G
INT5FK65	5F / 650 mm	grey	5F	18G
INT5FK70	5F / 700 mm	grey	5F	18G
INT5FK80	5F / 800 mm	grey	5F	18G
INT5FK90	5F / 900 mm	grey	5F	18G
INT6FK50	6F / 500 mm	green	6F	18G
INT6FK55	6F / 550 mm	green	6F	18G
INT6FK60	6F / 600 mm	green	6F	18G
INT6FK65	6F / 650 mm	green	6F	18G
INT6FK70	6F / 700 mm	green	6F	18G
INT6FK80	6F / 800 mm	green	6F	18G
INT6FK90	6F / 900 mm	green	6F	18G
INT7FK50	7F / 500 mm	orange	7F	18G
INT7FK55	7F / 550 mm	orange	7F	18G
INT7FK60	7F / 600 mm	orange	7F	18G
INT7FK65	7F / 650 mm	orange	7F	18G
INT7FK70	7F / 700 mm	orange	7F	18G
INT7FK80	7F / 800 mm	orange	7F	18G
INT7FK90	7F / 900 mm	orange	7F	18G
INT8FK50	8F / 500 mm	blue	8F	18G
INT8FK55	8F / 550 mm	blue	8F	18G
INT8FK60	8F / 600 mm	blue	8F	18G
INT8FK65	8F / 650 mm	blue	8F	18G
INT8FK70	8F / 700 mm	blue	8F	18G
INT8FK80	8F / 800 mm	blue	8F	18G
INT8FK90	8F / 900 mm	blue	8F	18G

Special order can contain introducer with radiopaque marker.
Model reference number: **INT8FK90M** – introducer with a marker.

ANAESTHESIOLOGICAL INTRODUCER

Set contains:

- sheath with a valve 8.5F x 11 cm
- angiographic needle 18G x 7 cm
- guide wire .035" / 40 cm / „J"
- dilatator 8.5F x 11 cm
- three-way stopcock
- sheath's protector
- syringe 10 ml

Catalog Number	Sheath size / length
INT8.5FA	5F / 110 mm

RADIAL ARTERY INTRODUCER FOR CARDIOLOGY

Set with hemostasis valve to introduce and exchange catheters and endocavitary electrodes

Set contains:

- sheath with a valve
- angiographic needle
- guide wire
- dilatator

Catalog Number	Sheath size / length	Cap colour	Dilatator size	Angiographic needle size	Guide wire size / length / type
INT5FKR	5F / 110 mm	grey	5F	21G	.018" / 45 cm / "J"
INT5FKRP	5F / 110 mm	grey	5F	21G	.018" / 45 cm / STRAIGHT
INT5F21KRP	5F / 110 mm	grey	5F	21G	.021" / 45 cm / STRAIGHT
INT5FKR23	5F / 230 mm	grey	5F	21G	.018" / 45 cm / "J"
INT5FKR23P	5F / 230 mm	grey	5F	21G	.018" / 45 cm / STRAIGHT
INT5F21KR23P	5F / 230 mm	grey	5F	21G	.021" / 45 cm / STRAIGHT
INT6FKR	6F / 110 mm	green	6F	21G	.018" / 45 cm / "J"
INT6FKRP	6F / 110 mm	green	6F	21G	.018" / 45 cm / STRAIGHT
INT6F21KRP	6F / 110 mm	green	6F	21G	.021" / 45 cm / STRAIGHT
INT6FKR23	6F / 230 mm	green	6F	21G	.018" / 45 cm / "J"
INT6FKR23P	6F / 230 mm	green	6F	21G	.018" / 45 cm / STRAIGHT
INT6F21KR23P	6F / 230 mm	green	6F	21G	.021" / 45 cm / STRAIGHT
INT7FKR	7F / 110 mm	orange	7F	21G	.018" / 45 cm / "J"
INT7FKRP	7F / 110 mm	orange	7F	21G	.018" / 45 cm / STRAIGHT
INT7F21KRP	7F / 110 mm	orange	7F	21G	.021" / 45 cm / STRAIGHT
INT7FKR23	7F / 230 mm	orange	7F	21G	.018" / 45 cm / "J"
INT7FKR23P	7F / 230 mm	orange	7F	21G	.018" / 45 cm / STRAIGHT
INT7F21KR23P	7F / 230 mm	orange	7F	21G	.021" / 45 cm / STRAIGHT

Special order can contain introducer with radiopaque marker.
 Model reference number: INT7FKRM – introducer with a marker.



The set contains a guide wire of special construction. The soft distal part of the guide wire enables its introduction to the vessel, and the stiff proximal part provides an excellent scaffold for the dilatator, along with preventing vessel perforation. The distal part of the sheath and dilatator is finished in a way that enables an easy and atraumatic introduction of the catheter to the vessel.

INTRODUCERS

BRAIDED INTRODUCER (COIL)

- High flexibility & trackability
- Kink-resistant
- Platinum coil marker – high visibility under fluoroscopy
- Specially designed atraumatic tip

Catalog Number	Sheath size / length	Dilatator size	Angiographic needle size	Use with guide wire
INT5F45T	5F / 450 mm	5F	18G	.035"
INT5F60T	5F / 600 mm	5F	18G	.035"
INT5F90T	5F / 900 mm	5F	18G	.035"
INT6F45T	6F / 450 mm	6F	18G	.035"
INT6F60T	6F / 600 mm	6F	18G	.035"
INT6F90T	6F / 900 mm	6F	18G	.035"
INT7F45T	7F / 450 mm	7F	18G	.035"
INT7F60T	7F / 600 mm	7F	18G	.035"
INT7F90T	7F / 900 mm	7F	18G	.035"
INT8F45T	8F / 450 mm	8F	18G	.035"
INT8F60T	8F / 600 mm	8F	18G </td <td>.035"</td>	.035"
INT8F90T	8F / 900 mm	8F	18G	.035"

Other diameters & other lengths (from 23 cm to 90 cm) are available on request.

Introducers are intended for introducing catheters and endocavitary electrodes.

The introducer's sheath is delivered with hemostasis valve, which prevents from losing blood during the surgery and it minimizes the embolism manifestation possibility.

The connection of the drain with sheath's head enables the infusion of fluids during the electrode placement procedure.

Caution:
This publication could include technical or other inaccuracies or typographical errors. Changes are periodically added to the information on herein; these changes will be incorporated in the new editions of the publication. Balton Ltd. may make improvements and/or changes in the product(s) described in this publication at any time. This publication is only a commercial information and cannot be handled the commercial offer in the meaning of the art. 66 §1 of Polish Civil Code.



Bioline Comercial Ltda.
R. José Getúlio, 579 14º andar cj 141/146
Liberdade | SP | 01509 001
+55 11 3274 7600
www.bioline.com.br